

# BPL\_YEAST\_AIR\_Fedbatch demo

This notebook demonstrate yeast fedbatch cultivation. We look at impact of changes in the glucose feeding. We also take a look at tuning of the DO-control system. Both liquid- and gasphase are included in the model.

Interaction with the compiled model as FMU is mainly through the simplified commands: `par()`, `init()`, `newplot()`, `simu()` etc. The last simulation is always available in the workspace and called 'sim\_res'. The command `describe()` brings mainly up description information from the actual Modelica code from the FMU but is complemented with information given in the dedicated Python setup-file.

The idea is to demonstrate how simulations and varying conditions can provide some process insight that can support the experimental work. I hope that at the end of this session you are ready to formulate your own questions you want to address with simulations - and you can just go on in this notebook! Just press the field "+Code" in the upper left part of notebook interface and you get a new "cell" where you write your own code. You can copy and paste from cells above using `ctrl-c` and `ctrl-p` as usual and edit the cell. When you are ready to execute the cell just press the "play button" to the left in the cell or press `shift-enter` as in "ordinary" Jupyter notebooks.

After a session you may want to save your own notebook. That you can do on your Google Drive account and I refer to Colab instructions for how to do this. It is easy.

Enjoy!

In [1]:

```
run -i BPL_YEAST_AIR_Fedbatch_DOcontrol_explore.py
```

Windows - run FMU pre-compiled JModelica 2.14

Model for bioreactor has been setup. Key commands:

- `par()` - change of parameters and initial values
- `init()` - change initial values only
- `simu()` - simulate and plot
- `newplot()` - make a new plot
- `show()` - show plot from previous simulation
- `disp()` - display parameters and initial values from the last simulation
- `describe()` - describe culture, broth, parameters, variables with values / units

Note that both `disp()` and `describe()` takes values from the last simulation

Brief information about a command by `help()`, eg `help(simu)`

Key system information is listed with the command `system_info()`

In [2]:

```
plt.rcParams['figure.figsize'] = [36/2.54, 30/2.54]
```

## About the process model

We can get information about the process, liquid- and gas-phase by the command `describe()`. This command can also be used to bring up information about a specific variable or parameter.

However, you should use `describe()` after a simulation to get the values used during the simulation.

```
In [3]: describe('culture'); print(); describe('liquidphase'); print(); describe('gasphase')
```

Saccharomyces cerevisiae - default parameters for strain H1022

Reactor broth substances included in the model

```
Cells index      = 1 - molecular weight = 24.6 Da
Glucose index    = 2 - molecular weight = 180.0 Da
Ethanol index    = 3 - molecular weight = 46.0 Da
Dissolved O2 index = 4 - molecular weight = 32.0 Da
Dissolved CO2 index = 5 - molecular weight = 44.0 Da
```

Reactor gasphase substances included in the model

```
N2 etc index = 1 - molecular weight = 28.0 Da
O2 index     = 2 - molecular weight = 32.0 Da
CO2 index    = 3 - molecular weight = 44.0 Da
Ethanol index = 4 - molecular weight = 46.0 Da
```

## First simulations - adjusting start of substrate feeding


```
In [4]: # Culture parameters and others at default values
        par(qO2lim=0.0069)

        # Process initial conditions
        init(V_0=4.5, VG_0=4.5*30, VX_0=4.5*2, VE_0=4.5*0)

        # Feed profile
        par(t_start=9, F_start=0.044, mu_feed=0.20, F_max=0.18)

        # DO-control parameters
        par(samplePeriod=1/60, K=10, Ti=0.5, I_0=500)

        # Simulate and plot
        newplot(title='Yeast fedbatch cultivation', plotType='Overview')
        simu(20)
```

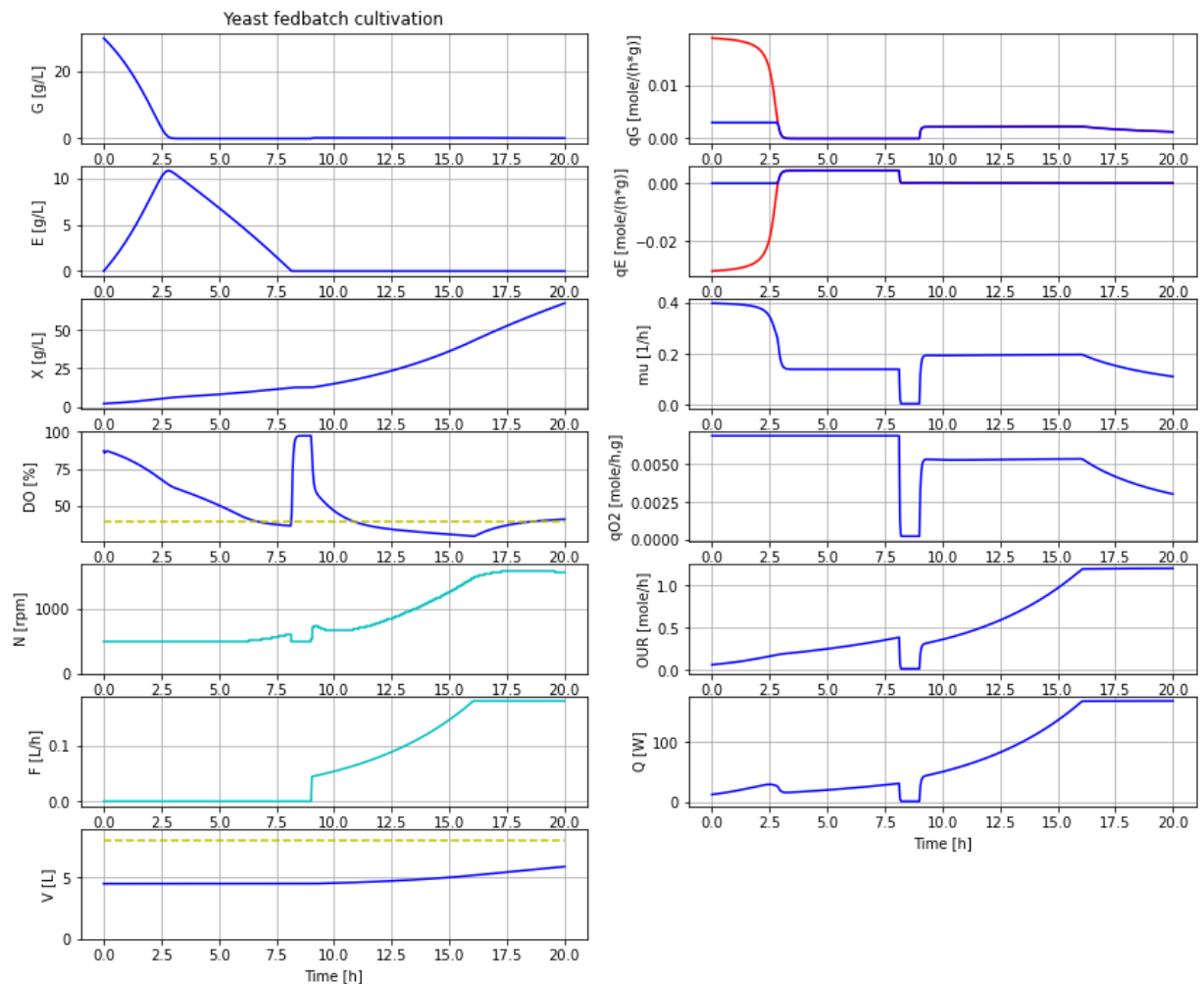

Now we can get value of broth volume as well as the headspace and values are the last ones in the simulation

```
In [5]: describe('bioreactor.V')
```

Reactor broth volume : 5.892 [ L ]

```
In [6]: describe('bioreactor.V_gasphase')
```


Volume of the gas phase : 2.108 [ L ]

```
In [7]: # Take a look at the parameters available to adjust the dosage scheme
disp('dosage', decimals=4)
```

```
mu_feed : 0.2
F_0 : 0.0
t_start : 9.0
F_start : 0.044
F_max : 0.18
```

```
In [8]: # Let us start the feeding just after the batch phase has ended and keep other param
par(t_start=8.1)

# Simulate and plot
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
simu(20)
```

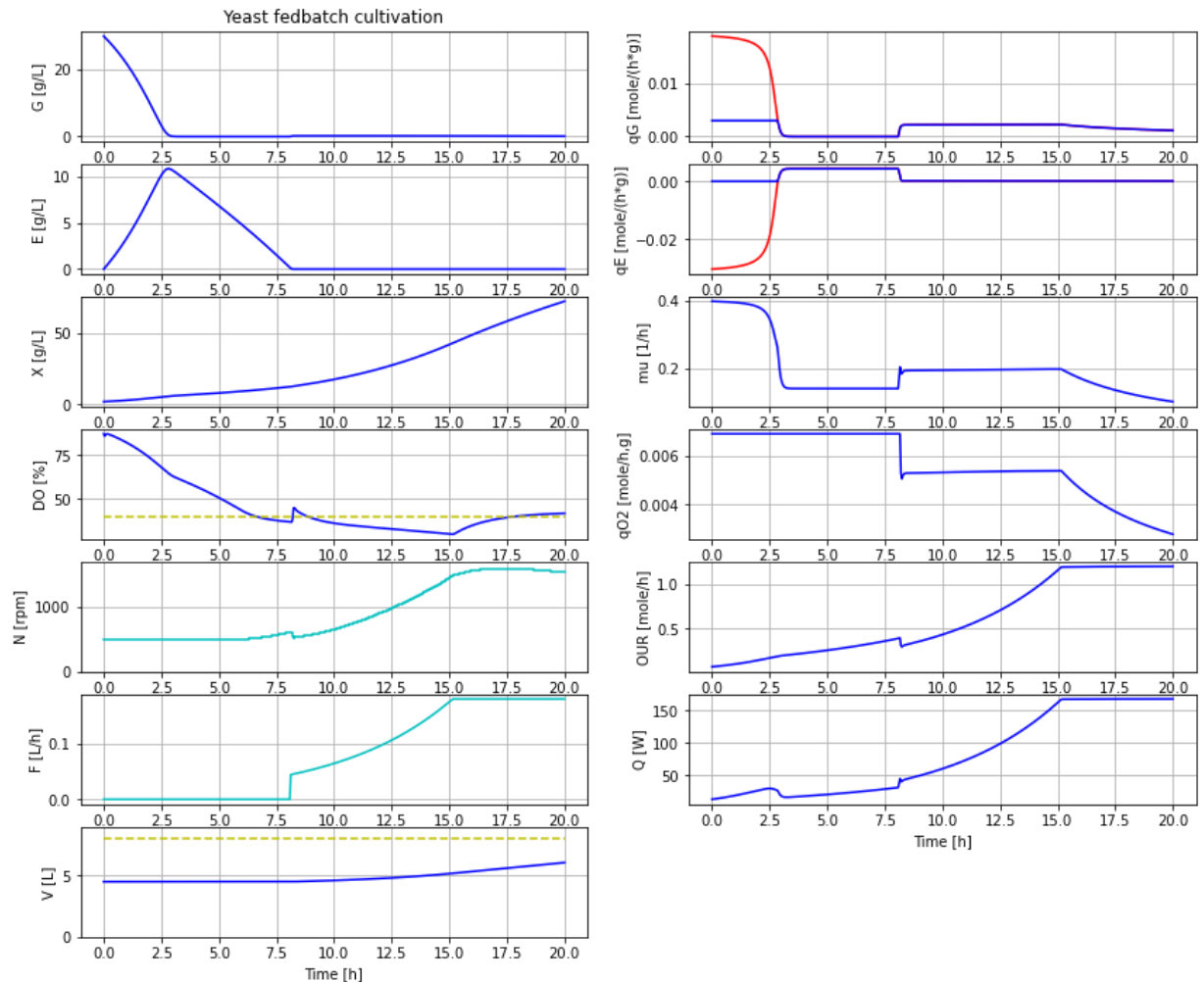


The increase of DO to about 50 % at end of batch phase should be possible to detect easily. This simulation is more realistic and we use these settings from now on.

## DO-control - tuning of PI-regulator parameters

Let us focus on the DO-control system and choose a more limited plotType. We study the impact of PI control parameters and see if we can decrease the control error without loosing stability.

```
In [9]: # Let us take a closer look at the DO-control system and try to make control error s
newplot(title='Yeast fedbatch cultivation - DO-control - increase K', plotType='Focu
for value in [10, 20, 40]: par(K=value); simu(20)

# Reset K to the original value
par(K=10)
```

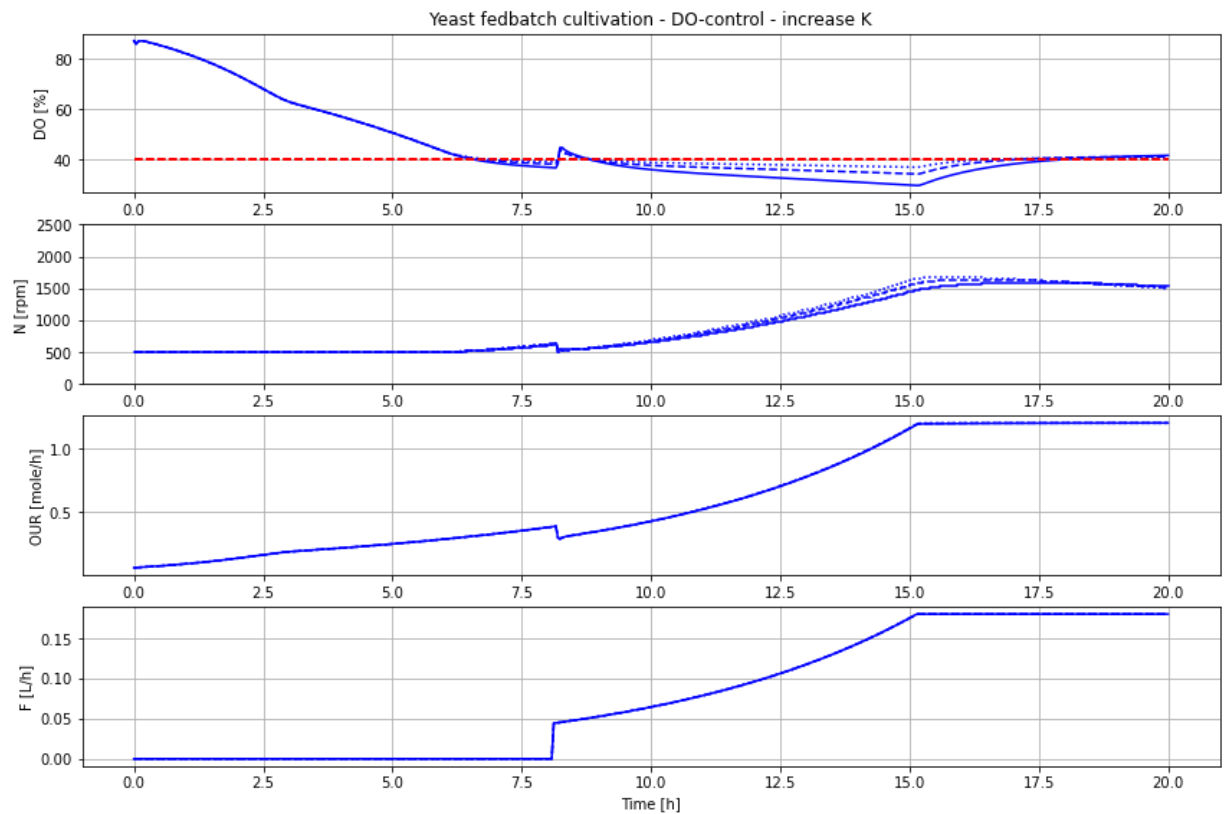


We see that by a higher control gain  $K$  the DO-control error get smaller and the stability of the control system is maintained.

**Exercise** I leave for you to study the impact variation of the  $T_i$ -parameter. Just make a new cell below. Then copy and paste the cell above and change parameter to  $T_i$ .

## Sensitivity to changes in feed-profile


Now, let us focus on investigating impact of changes in the feed-profile. The goal is to increase the produced cell mass without accumulation of by-product ethanol. Simulation can bring some insight into how behaviour of the different variables change when by-product is formed. This insight can help to interpret experimental results.

```
In [10]: # Let us check the sensitivity to changes in the feed profile design
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
for value in [0.044, 0.038, 0.050]: par(F_start=value); simu(20)
```

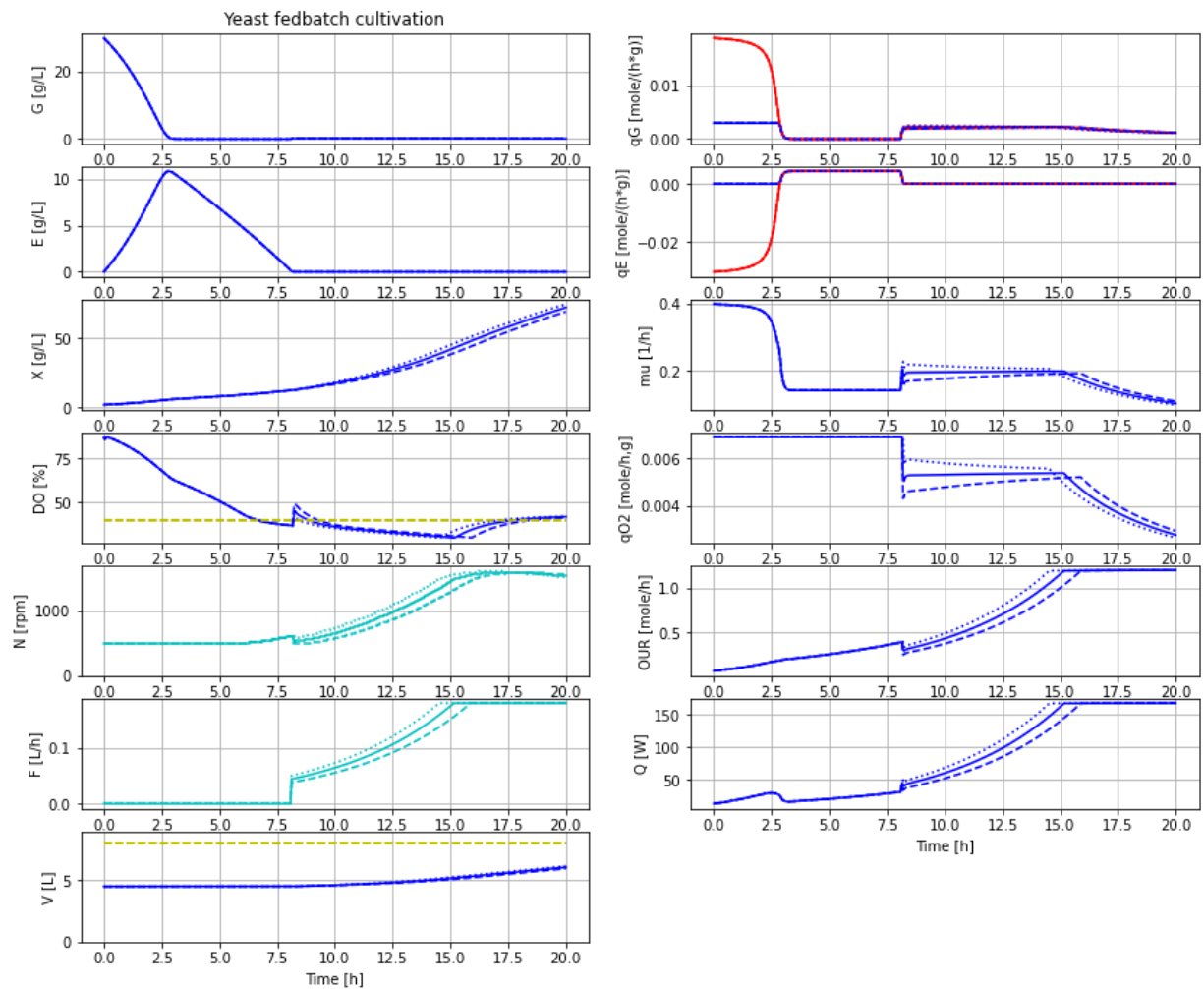


The variation in  $F_{\text{start}}$  has an impact and we see that the actual growth rate during fedbatch phase do converge to the set growth rate of the feed, but it takes more than 5 hours.

```
In [11]: # Let us investigate a feedprofile that is closer to the maximal capacity
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
par(F_start=0.044, mu_feed=0.20); simu(20)
par(F_start=0.050, mu_feed=0.22); simu(20)
par(F_start=0.057, mu_feed=0.26); simu(20)
```

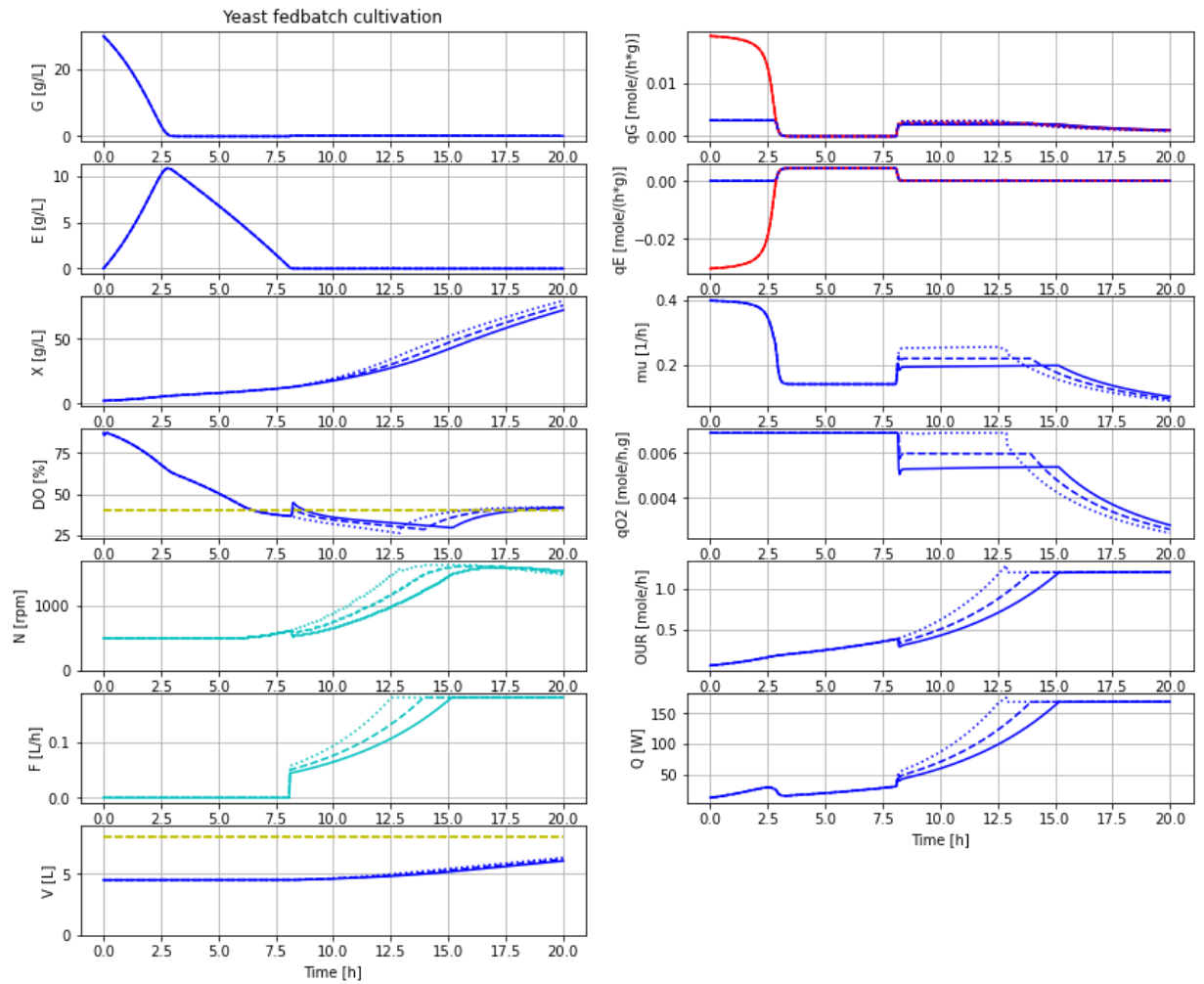

In [12]:

```
# And let us see what happens if the feedprofile exceed the culture capacity
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
par(F_start=0.057, mu_feed=0.26); simu(20)
par(F_start=0.063, mu_feed=0.28); simu(20)
par(F_start=0.044, mu_feed=0.20)
```

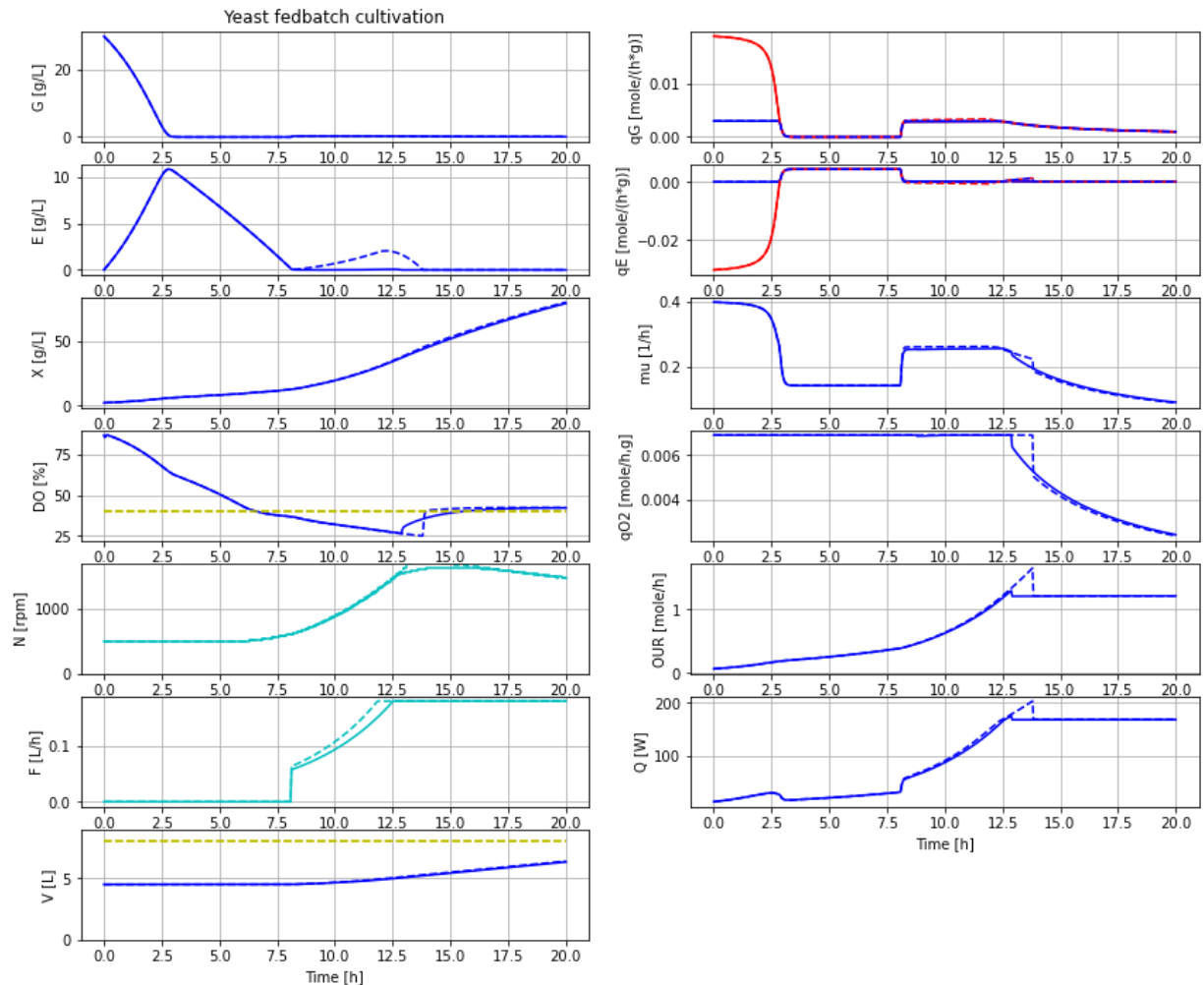

Note that with the feedprofile that exceed culture respiratory capacity, ethanol is accumulated during time 8-12.5 hours. When the feedprofile then is constant from time 12.5 hours and on, then the accumulated ethanol is consumed over about an hour. This leads to a higher oxygen demand and heat production during this time. The specific cell growth rate is also slightly higher during this period.

**Exercise** You can investigate the impact of changing the maximal feedrate  $F_{\max}$ . Make sure that the DO level do not get too low.

## Make your own diagrams

There are a couple of pre-defined plotType for the application that you make by the command `newplot()`. The command result in a list "diagrams" that describe the commands that make the plot when you call `simu()` or you just want to look at the last simulation again with a changed plotType using `show()`.

You can also in Jupyter notebook directly define the list "diagrams" and then that will be used for subsequent calls of `simu()` or `show()`. When you have made a diagram that you want to reuse many times you can bring it into the python-setup file and edit the `newplot()` command and add a new plotType.

Below a few simple examples that show how to do a diagram directly in the notebook

```
In [13]: # Improvise and make your own diagram - cell concentration in a logarithmic plot
plt.figure()
```




```

ax1 = plt.subplot(1,1,1)
ax1.set_ylabel('X [g/L]')
ax1.set_xlabel('Time [h]')
ax1.grid()

setLines()
diagrams.clear()
diagrams.append("ax1.semilogy(sim_res['time'], sim_res['bioreactor.c[1]'], color='b'
show()

```

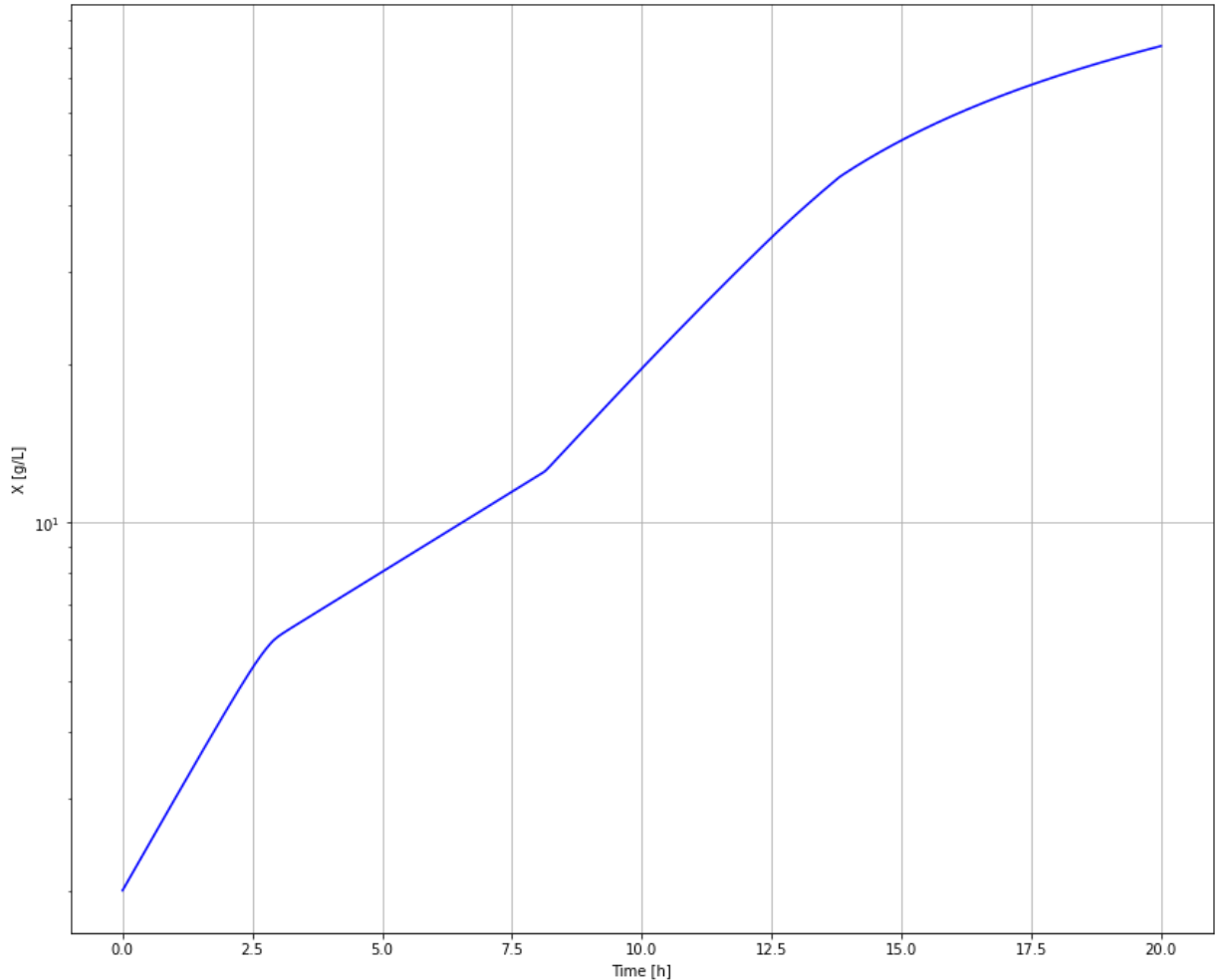


In [14]:

```

# - study the variation of Kla together with DO and N during cultivation
# - study the relation qO2 vs qG(G)
plt.figure()
ax1 = plt.subplot(3,1,1); ax2 = plt.subplot(3,1,2); ax3 = plt.subplot(3,1,3)
ax1.set_ylabel('DO [%]'); ax1.grid()
ax2.set_ylabel('Kla [1/h]'); ax2.grid()
ax3.set_ylabel('N [rpm]'); ax3.grid()
ax3.set_xlabel('Time [h]')

setLines()
diagrams.clear()
diagrams.append("ax1.plot(sim_res['time'], sim_res['DOsensor.out'], color='b', lines
diagrams.append("ax2.plot(sim_res['time'], sim_res['bioreactor.gas_liquid_transfer.K
diagrams.append("ax3.plot(sim_res['time'], sim_res['bioreactor.N'], color='c', lines
show()

```

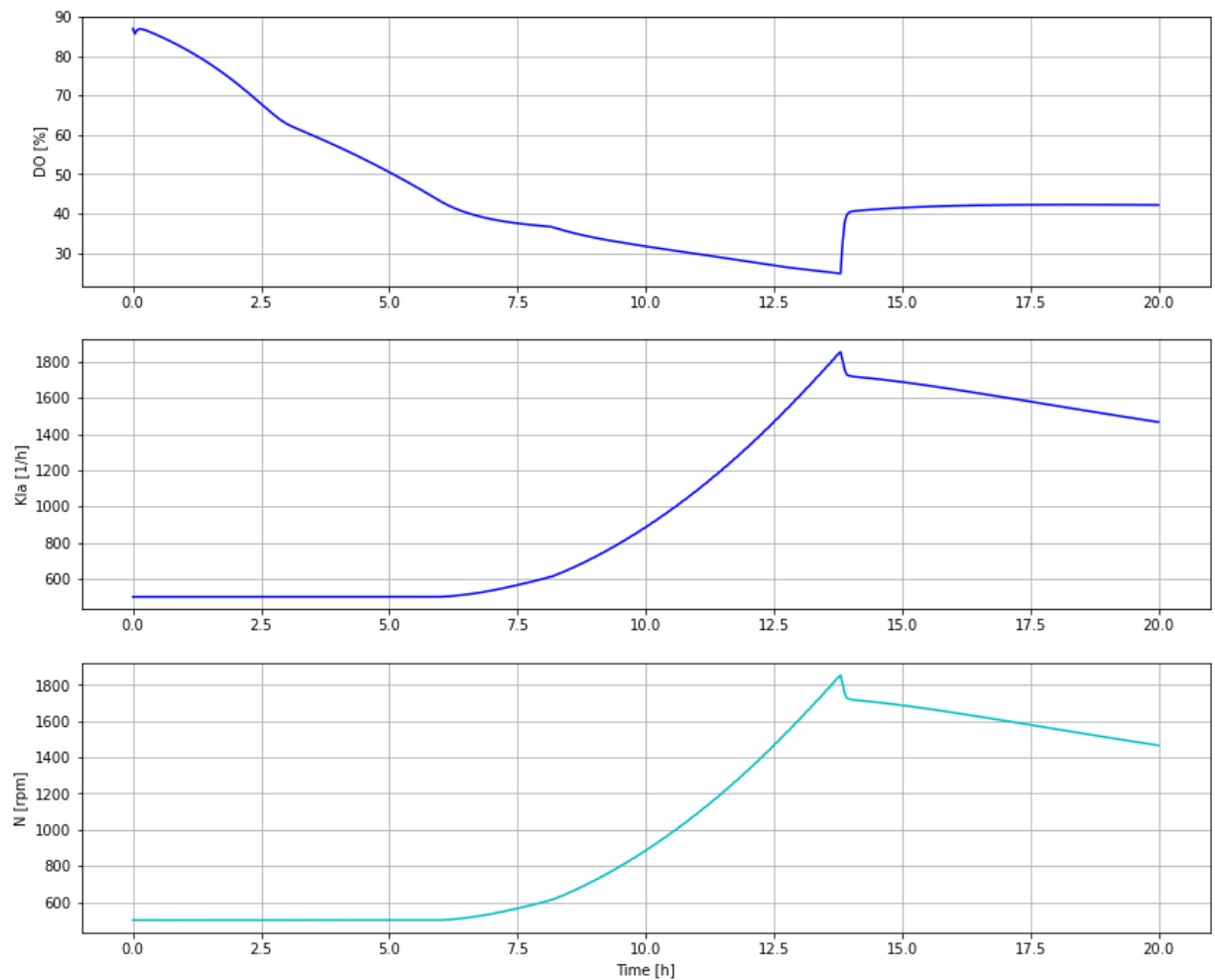


The relation is  $KLa = \alpha_{O2} \cdot N$  and we see the value of the parameter should be around 1.0, and we check below

```
In [15]: disp('bioreactor.gas_liquid_transfer.alpha_O2')
```

```
alpha_O2 : 1.0
```

```
In [16]: # - study the relation qO2 vs qG(G)
plt.figure()
ax1 = plt.subplot(1,1,1)
ax1.set_ylabel('qO2 [ ]')
ax1.set_xlabel('qG [ ]')
ax1.grid()

setLines()
diagrams.clear()
diagrams.append("ax1.plot(sim_res['bioreactor.culture.qGm'], sim_res['bioreactor.cul
par(F_start=0.057, mu_feed=0.26)
simu(20)
```

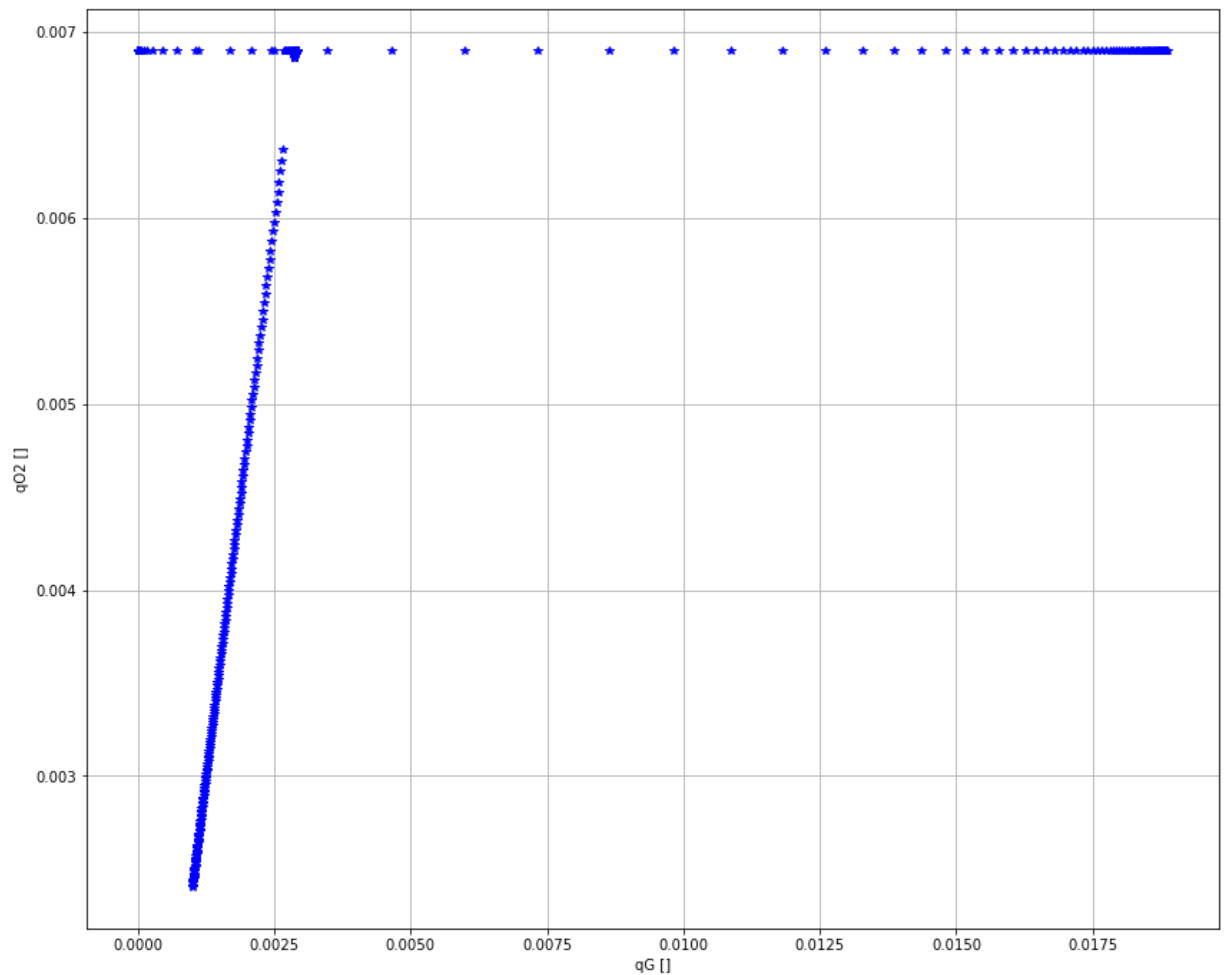

During the cultivation we have a number of data points for qG and qO2 at the same time, during different conditions. What we see in the diagram is that qO2 increase with qG until qG reach a level of just above 0.0025 and then qO2 saturats for highter qG. This what expect to see.

We also see that for lower qG we have also qO2 values at saturation level. This points correspond to a situation where ethanol is consumed with the remaining respiratory capacity. Glucose is consumed by priority.

```
In [17]: model.get('bioreactor.outlet_gas.x_gas[4]')
```


```
Out[17]: array([-0.00713861])
```

```
In [18]: # Improvise and make your own diagram - study gasphase
plt.figure()
ax1 = plt.subplot(1,1,1)
ax1.set_ylabel('x_gas[02] []')
ax1.set_xlabel('Time [h]')
ax1.grid()

setLines()

diagrams.clear()
diagrams.append("ax1.plot(sim_res['time'], sim_res['bioreactor.outlet_gas.x_gas[4]'])
#diagrams.append("ax1.set_ylim([-200,0])")

show()
```

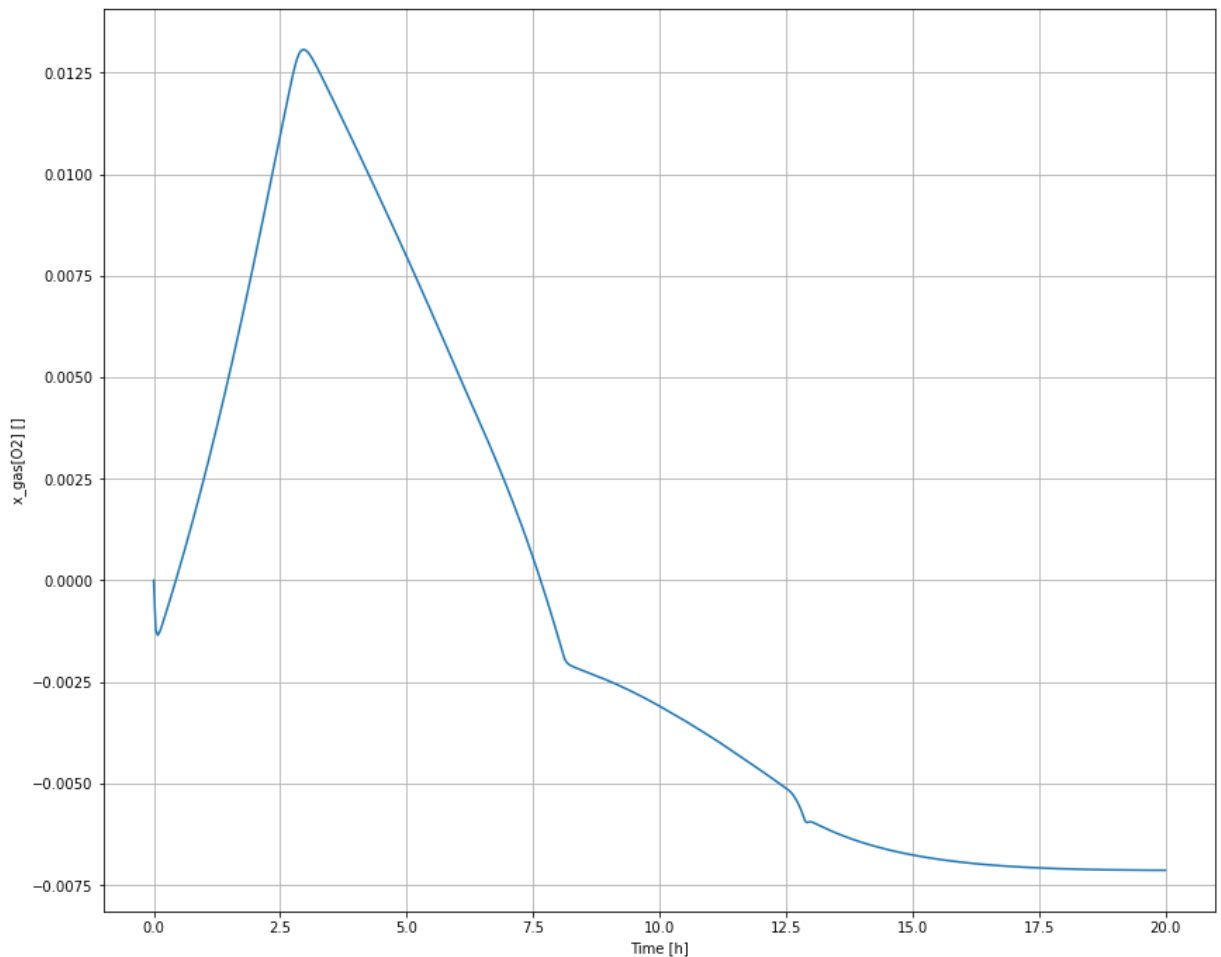

## Summary

- We have first seen an overview diagram of a typical yeast fedbatch cultivation where the feed started about an hour after the batch phase was finished. A new simulation was made where the feed started directly after detection of lack of substrate.
- We also took a look at the DO-control system and saw that we could decrease the control error by increasing the PI-controller gain. Stability of the control system remained.
- Then we tested variations in the feed dosage scheme and investigated the possibilities to increase the production.
- We also saw what happens if the feed dosage exceed the culture respiratory capacity and what to look for during the experimental work.
- Finally we saw some examples of how to improvise new diagrams.

## Appendix

```
In [19]: # List of components in the process setup and also a couple of other things like Liq
describe('parts')
```

```
['airFlow_setpoint', 'airtube', 'atmosphere', 'bioreactor', 'bioreactor.culture', 'bioreactor.gas_liquid_transfer', 'compressor', 'DO_setpoint', 'dosagescheme', 'DOsensor', 'feedtank', 'gasphase', 'liquidphase', 'MSL', 'N_high', 'N_low', 'PIreg', 'pump']
```

```
In [20]: system_info()
```

System information

- OS: Windows
- Python: 3.9.5
- PyFMI: 2.9.5
- FMU by: JModelica.org
- FMI: 2.0
- Type: FMUModelCS2
- Name: BPL\_YEAST\_AIR.Fedbatch\_DOcontrol
- Generated: 2022-08-19T11:13:22
- MSL: 3.2.2 build 3
- Description: Bioprocess Library version 2.1.0 beta
- Interaction: FMU-explore ver 0.9.2

In [ ]: